

Product Specifications




CE₀₁₂₃

Type	Specifications	Size (Fr)	Effective Length (mm)	Specifications	Size (Fr)	Effective Length (mm)	Specifications	Size (Fr)	Effective Length (mm)	Specifications	Size (Fr)	Effective Length (mm)
B	QB040450A	4	450	QB060900A	6	900	QB100450A	10	450	QB060900C	6	900
	QB040550A	4	550	QB061100A	6	1100	QB120450A	12	450	QB070450C	7	450
	QB040700A	4	700	QB070450A	7	450	QB060400B	6	400	QB070550C	7	550
	QB040900A	4	900	QB070550A	7	550	QB070400B	7	400	QB070700C	7	700
	QB041100A	4	1100	QB070700A	7	700	QB080400B	8	400	QB080450C	8	450
	QB050450A	5	450	QB070800A	7	800	QB040550C	4	550	QB080550C	8	550
	QB050550A	5	550	QB070900A	7	900	QB040700C	4	700	QB080700C	8	700
	QB050700A	5	700	QB071100A	7	1100	QB050450C	5	450	QB090450C	9	450
	QB050900A	5	900	QB080450A	8	450	QB050550C	5	550	QB100450C	10	450
	QB051100A	5	1100	QB080550A	8	550	QB050700C	5	700	QB120450C	12	450
	QB060450A	6	450	QB080700A	8	700	QB050900C	5	900	QB050450D	5	450
	QB060550A	6	550	QB090450A	9	450	QB060450C	6	450	QB060450D	6	450
	QB060700A	6	700	QB090550A	9	550	QB060550C	6	550	QB070450D	7	450
	QB060800A	6	800	QB090700A	9	700	QB060700C	6	700	—	—	—
	QC040450A	4	450	QC060900A	6	900	QC100450A	10	450	QC060900C	6	900
	QC040550A	4	550	QC061100A	6	1100	QC120450A	12	450	QC070450C	7	450
C	QC040700A	4	700	QC070450A	7	450	QC060400B	6	400	QC070550C	7	550
	QC040900A	4	900	QC070550A	7	550	QC070400B	7	400	QC070700C	7	700
	QC041100A	4	1100	QC070700A	7	700	QC080400B	8	400	QC080450C	8	450
	QC050450A	5	450	QC070800A	7	800	QC040550C	4	550	QC080550C	8	550
	QC050550A	5	550	QC070900A	7	900	QC040700C	4	700	QC080700C	8	700
	QC050700A	5	700	QC071100A	7	1100	QC050450C	5	450	QC090450C	9	450
	QC050900A	5	900	QC080450A	8	450	QC050550C	5	550	QC100450C	10	450
	QC051100A	5	1100	QC080550A	8	550	QC050700C	5	700	QC120450C	12	450
	QC060450A	6	450	QC080700A	8	700	QC050900C	5	900	QC050450D	5	450
	QC060550A	6	550	QC090450A	9	450	QC060450C	6	450	QC060450D	6	450
	QC060700A	6	700	QC090550A	9	550	QC060550C	6	550	QC070450D	7	450
	QC060800A	6	800	QC090700A	9	700	QC060700C	6	700	—	—	—
	QD040900A	4	900	QD060900A	6	900	QD080900A	8	900	QD070550C	7	550
	QD041100A	4	1100	QD061100A	6	1100	QD060450B	6	450	QD050550D	5	550
	QD050900A	5	900	QD070800A	7	800	QD070450B	7	450	QD060550D	6	550
	QD051100A	5	1100	QD070900A	7	900	QD050550C	5	550	QD070550D	7	550
	QD060800A	6	800	QD071100A	7	1100	QD060550C	6	550	—	—	—

HEETY™

Catheter Sheath Sets

Type B/C/D

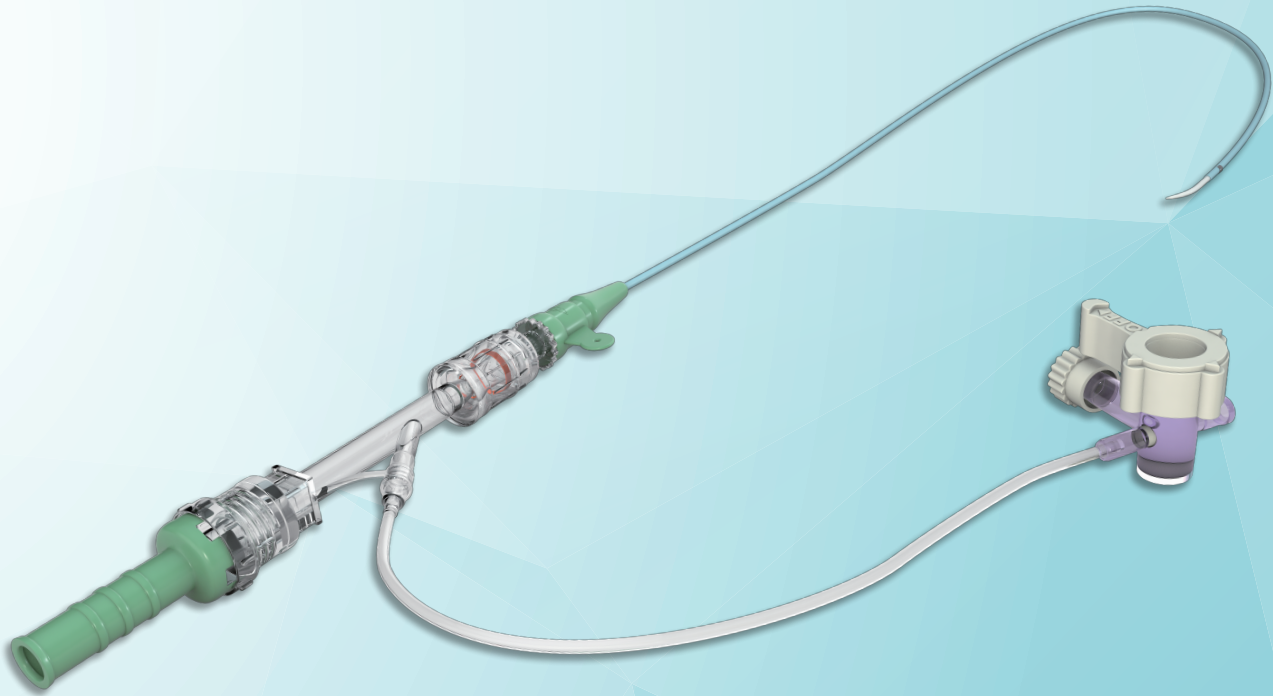

Zhejiang Barty Medical Technology Co.,Ltd.


Add : South 2F,Building 2-1,No.20 Street,ETDZ, Hangzhou,Zhejiang,China
 Tel : 0571-883737191
 Fax : 0571-88373716
 Web: //www.bartymedical.com www.vascularbarty.com
 Email: info@bartymedical.com

HEETY™

Catheter Sheath Sets

Type B/C/D

Variety Of Specifications And Types For More Applications

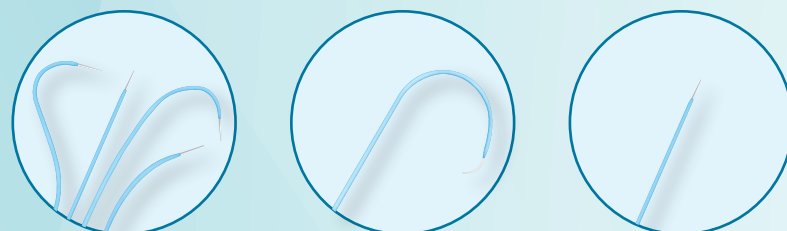


Type D: Y-connector

Design of Y-valve prevents aeroembolism, reduces blood leak, being convenient for operation. No Y-connector in Type A/B/C.


Excellent Flexural-resistance And compliance

Excellent flexural-resistance and compliance:
PTFE inner layer ensures smooth advancement;
Stainless steel inner layer ensures excellent flexural-resistance;
PA outer layer ensures smooth advancement inside blood vessel.

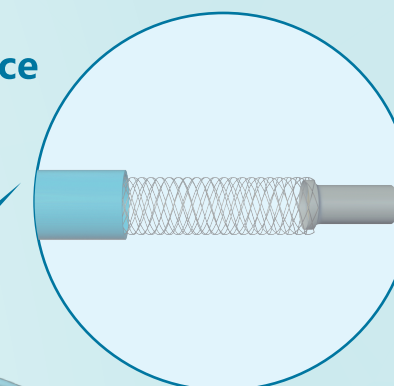

Hydrophilic Coating For Easy Insertion

Hydrophilic coating over sheath surface to provide a smooth surface for insertion.

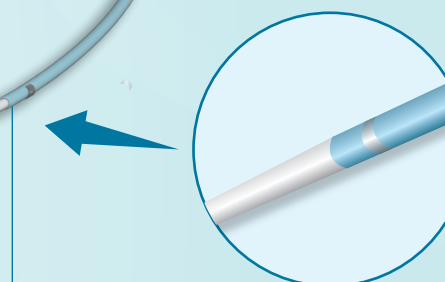

Safe And Humanized Ddesign To Lower Operational Risk

Perfect matching between sheath and dilator to reduce damage to blood vessel during operation;
Mark band with non transparent X ray on sheath distal end to locate.

Type	Product Configuration					
B	Catheter Sheath	Dilator	—	—	—	—
C	Catheter Sheath	Dilator	Introducer Needle	Scalpel blade	Guide Wire	Syringe
D	Catheter Sheath (Y- connector)	Dilator	—	—	—	—

Application

This product is intended to percutaneously insert into blood vessel during interventional operation, providing channel for angiographic catheter, guide catheter, balloon dilatation catheter and other similar devices into artery and vein.